



BLACK  
FOREST  
MEDICAL  
GROUP

# DORO LUCENT® IMRI HEADREST SYSTEMS



# DORO LUCENT® iMRI Headrest Systems



## FAST

Directly mountable to transfer boards of selected OR table manufacturers.

## FLEXIBLE

The parallelogram design supports you in properly positioning your patient for surgery.

## FAMILIAR

The classic 3-pin skull clamp design with ratcheting system allows you to rely on your trained routine and known workflows.

## AUTOMATIC IMAGE REGISTRATION (AIR) READY

Brainlab Automatic Image Registration (AIR) is the fundamental technology for high-precision surgical navigation from start to finish. It adds intraoperatively acquired MR images to surgical navigation while maintaining the preoperative surgical plan. This updated registration compensates for potential loss of navigation accuracy during surgery. The Brainlab registration matrix with integrated MR markers is designed for use exclusively with DORO LUCENT® iMRI Headrest System and therefore does not compromise standard patient positioning.





Please find further information on the product's web page.

# Advanced Technologies, Known Routines

Intra-operative Magnetic Resonance Imaging (iMRI) helps surgeons capture views of the brain during surgery. This allows them to remove a tumor more safely having the visual proof that it is completely eliminated. Our aim is to improve surgical procedures. Therefore, we developed a technically advanced solution, which can be used for iMRI, but is yet easy to use.


## DORO LUCENT® HEADPLATE INCLUDED



The DORO LUCENT® Headplate optimizes the clinical workflow of the DORO LUCENT® iMRI Headrest Systems as it grants the clinical staff the opportunity to position the patient as required for surgery before starting the actual pinning procedure, rendering later on shifts of the sedated patient obsolete. This is less procedural risk, less stress for the patient and more convenience for the surgeon and the OR team.


## MOVEMENT CAPABILITIES

DORO LUCENT® Parallelogram Adaptors allow you to shift, tilt, angle and rotate the skull clamp to always find a fitting position for your patient.

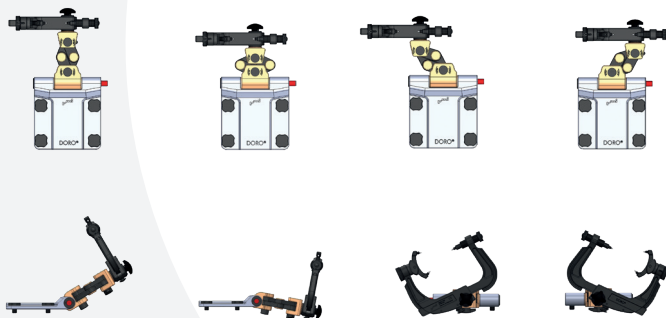

# DORO LUCENT® iMRI HEADREST SYSTEMS

- MRI SCANNER / OR TABLE COMBINATION OVERVIEW -

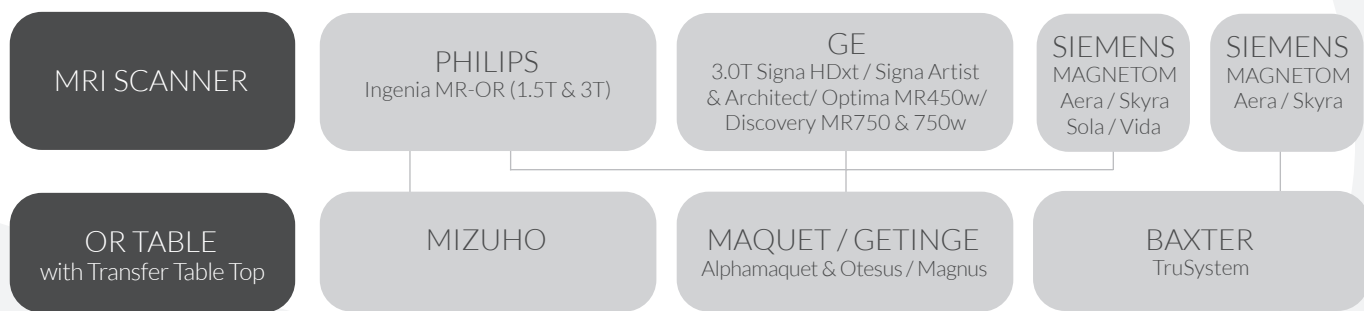

DORO LUCENT® iMRI Headrest System MIZUHO Item no. 4003.350		DORO LUCENT® iMRI Headrest System MAQUET Item no. 4003.300		DORO LUCENT® iMRI Headrest System Baxter Item no. 4003.200	
Included Items	Item no.	Included Items	Item no.	Included Items	Item no.
DORO® Skull Clamp Radiolucent for DORO LUCENT® Parallelogram Adaptor	4003.001	DORO® Skull Clamp Radiolucent for DORO LUCENT® Parallelogram Adaptor	4003.001	DORO® Skull Clamp Radiolucent for DORO LUCENT® Parallelogram Adaptor	4003.001
DORO® Headplate	4003.005	DORO® Headplate	4003.005	DORO® Headplate	4003.005
DORO® Adjusting Wrench for DORO LUCENT® Parallelogram Adaptor	4003.011	DORO® Adjusting Wrench for DORO LUCENT® Parallelogram Adaptor	4003.011	DORO® Adjusting Wrench for DORO LUCENT® Parallelogram Adaptor	4003.011
DORO LUCENT® Storage Case for Parallelogram Adaptors	4003.030	DORO LUCENT® Storage Case for Parallelogram Adaptors	4003.030	DORO LUCENT® Storage Case for Parallelogram Adaptors	4003.030
DORO LUCENT® iMRI Parallelogram Adaptor MAQUET	4003.311	DORO LUCENT® iMRI Parallelogram Adaptor MAQUET	4003.311	DORO LUCENT® iMRI Parallelogram Adaptor Baxter	4003.211
DORO LUCENT® Disposable Skull Pins, Adult (3 pcs)	1106.003	DORO LUCENT® iMRI Transfer Collision Indicator	4003.XXX	DORO LUCENT® iMRI Pad Baxter	4003.212
A fitting iMRI Transfer Collision Indicator must be purchased from Phillips.		DORO LUCENT® Disposable Skull Pins, Adult (3 pcs)	1106.003	DORO LUCENT® iMRI Transfer Collision Indicator	4003.XXX
				DORO LUCENT® Disposable Skull Pins, Adult (3 pcs)	1106.003
				DORO® Torx Screwdriver T45, nonmagnetic	3033-68

**PLEASE NOTE:**

With the first installation of a DORO LUCENT® iMRI Headrest System, an expert training is mandatory. Please contact your local DORO® representative for details.



DORO LUCENT® Headplate  
Item no. 4003.005



DORO LUCENT® Disposable Skull Pin  
Item no. 1106.003



DORO® Adjusting Wrench for Parallelogram Adaptor  
Item no. 4003.011



DORO LUCENT® iMRI Transfer Collision Indicator  
Item no. 4003.XXX



DORO LUCENT® Storage Case for Parallelogram Adaptors  
Item no. 4003.030



pro med instruments – part of Black Forest Medical Group  
Headquarters | Bötzingen Straße 86 | 79111 Freiburg, Germany  
Phone +49 761 384 222 10 | Email info@blackforestmedical.com

US Subsidiary | 4529 S.E. 16th Place, Suite 101 | Cape Coral, FL 33904  
Toll Free 877 225 4086 | Email info.us@blackforestmedical.com

Version 07 © 2024 pro med instruments GmbH, Freiburg, Germany.

All rights reserved. The information provided was compiled with utmost care. Errors can, however, not be excluded. Errors do not entitle to any legal claims. We reserve the right to modify the products at any time.