

Customer: _____

Contact Person: _____

Address: _____

Telephone: _____

Email: _____



Date:

TEST REPORT – COMPATIBILITY WITH DORO® SKULL PINS

The following skull clamps were tested by the undersigned for compatibility with DORO® SKULL PINS, manufactured by pro med instruments, part of Black Forest Medical Group (herein after BFMG).

Test Procedure

The test gauge has a “go” and a “no go” side (red indication shows the “no go” side). Each receptacle bore of the skull clamp must be tested with both sides. Before testing, please ensure that the receptacles are clean and not damaged. For such test, the test gauge must be introduced without any force (only by its own weight):

- The “No Go” side must not fit into the receptacle
- The “Go” side shall move into the receptacle without any additional force

✓ **If both criteria are fulfilled**, compatibility of DORO® SKULL PINS to the reviewed third-party product under consideration is ensured.

✗ **If one of these 2 criteria is not met**, the reviewed third-party product is not compatible to the DORO® SKULL PINS.


Results shall be documented below and the BFMG representative shall archive this document.



DORO® Skull Pin types



DORO® Pin Receptacles

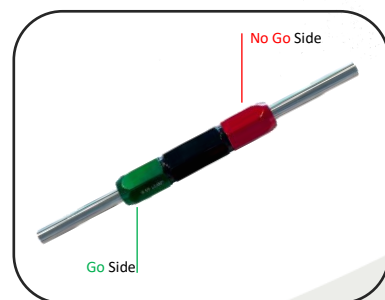

DORO® Skull Pin Limit Gauge

Customer: _____

Contact Person: _____

Address: _____

Telephone: _____

Email: _____



BLACK
FOREST
MEDICAL
GROUP

Products Under Consideration

Manufacturer					
Item number					
Serial number					
Acquisition date					
2- pin receptacle	<input type="checkbox"/> GO <input type="checkbox"/> NO GO <input type="checkbox"/> GO <input type="checkbox"/> NO GO	<input type="checkbox"/> GO <input type="checkbox"/> NO GO <input type="checkbox"/> GO <input type="checkbox"/> NO GO	<input type="checkbox"/> GO <input type="checkbox"/> NO GO <input type="checkbox"/> GO <input type="checkbox"/> NO GO	<input type="checkbox"/> GO <input type="checkbox"/> NO GO <input type="checkbox"/> GO <input type="checkbox"/> NO GO	<input type="checkbox"/> GO <input type="checkbox"/> NO GO <input type="checkbox"/> GO <input type="checkbox"/> NO GO
1- pin receptacle	<input type="checkbox"/> GO <input type="checkbox"/> NO GO	<input type="checkbox"/> GO <input type="checkbox"/> NO GO	<input type="checkbox"/> GO <input type="checkbox"/> NO GO	<input type="checkbox"/> GO <input type="checkbox"/> NO GO	<input type="checkbox"/> GO <input type="checkbox"/> NO GO
Comments					
Compatibility with DORO® SKULL PINS	<input type="checkbox"/> YES <input type="checkbox"/> NO	<input type="checkbox"/> YES <input type="checkbox"/> NO	<input type="checkbox"/> YES <input type="checkbox"/> NO	<input type="checkbox"/> YES <input type="checkbox"/> NO	<input type="checkbox"/> YES <input type="checkbox"/> NO

Test Equipment Name	Number	Valid Until

Signature BFMG representative / date

Signature customer / date

pro med instruments GmbH
part of Black Forest Medical Group
Bötzingler Straße 86
79111 Freiburg im Breisgau, Germany
Telefon: +49 761 384 222 10
Fax: +49 761 384 222 80
Email: info@blackforestmedical.com

Volksbank Freiburg eG
IBAN: DE20 6809 0000 0051 7702 18
SWIFT-BIC: GENO DE 61 FR1

Commerzbank Freiburg
IBAN: DE44 6808 0030 0402 0182 00
SWIFT-BIC: DRES DE FF 680

HypoVereinsbank Freiburg
IBAN: DE48 6802 0186 0321 3114 58
SWIFT-BIC: HYVE DE MM 357

Geschäftsführer
Matthias E. Schüle

Registergericht Freiburg
HRB 48 79

Ust. Nr. / VAT Reg. No.
DE 811 778 789

www.blackforestmedical.com